



810 LX Washer

Undercounter or freestanding

- 1 to 2 levels
- Upper and lower spray arm
- Steam condenser

Dimensions:

Chamber: 20"H x 21"W x 20"D

External: 33.25"H x 23.5"W x 26.75"D



910 LX Washer/Dryer

Freestanding

- HEPA filtered chamber and direct injector drying
- 1 to 2 level
- Upper and lower spray arm

Dimensions:

Chamber: 21.40"H x 21"W x 20"D

External: 57.75"H x 23.75"W x 29"D



1600 UP Washer/Dryer

Freestanding

- Chamber drying standard (injector drying on SS models)
- 1 to 4 levels
- Upper and 5 lower spray arms

Dimensions:

Chamber: 33.75"H x 25.75"W x 26"D

External: 70"H x 34.75"W x 33.75"D



815 LX Washer/Dryer

Undercounter or freestanding

- Forced air drying system
- 1 to 2 levels
- Upper and lower spray arm
- Steam condenser

Dimensions:

Chamber: 20"H x 21"W x 20"D

External: 33.25"H x 23.5"W x 26.75"D



1300 LX Washer/Dryer

Freestanding

- HEPA filtered chamber and direct injector drying
- 1 to 3 levels
- Upper and lower spray arm

Dimensions:

Chamber: 25.35"H x 21"W x 20"D

External: 64"H x 23.75"W x 29"D



1800 LXA Washer/Dryer

Freestanding

- HEPA filtered chamber and direct injector drying
- 1 to 6 levels
- Automatic vertical, sliding glass door
- Upper and lower spray arms

Dimensions:

Chamber: 30.3"H x 41.9"W x 20"D

External: 78"H x 52.6"W x 32.5"D



820 LX Washer/Dryer

Undercounter or freestanding

- HEPA filtered chamber and direct injector drying
- 1 to 2 levels
- Upper and lower spray arm
- Steam condenser

Dimensions:

Chamber: 20"H x 21"W x 20"D

External: 33.25"H x 35.5"W x 26.75"D



1400 LXP Washer/Dryer

Freestanding

- HEPA filtered chamber and direct injector drying
- 1 to 3 levels
- Upper and lower spray arm

Dimensions:

Chamber: 30.3"H x 24.4" W x 24.4"D

External: 68"H x 31"W x 32.7"D



TI Series Washer/Dryers

Freestanding

- Industrial applications
- Multiple chamber sizes
- Single door or pass through design

Contact LANCER for further details on specific models.

1600 TISS shown

Racks and Accessories

Basic Baskets

Used for cleaning beakers and other large necked items on lowest level.

Basic basket with an added spray arm is available for cleaning on upper levels.



PST - Basket

Injector Racks

For cleaning and drying inside narrow necked laboratory glassware using injectors. Spindle headers inject water, detergent and acid directly into glassware.



IXL - Jet Rack

Jet racks are available in various combinations:


- 24 to 144 jets (based on washer model)
- Long, short or combination of jets
- Various diameters - (2mm, 6mm, 8mm, combination)

Other Racks/Accessories

- Pipette Racks (grommet and canister styles)
- DNA Sequencing Plate Rack
- Culture Dish Rack
- Mesh Baskets (for cleaning small items)
- Combination Racks
- Loading Trolleys/Carts
- Stainless Steel Base Cabinets (undecounter models)



PPI - Pipette Rack

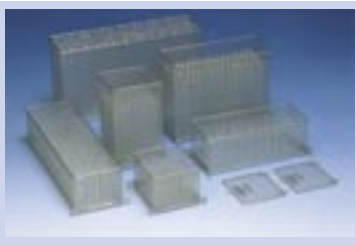

Mesh Baskets

Rack customization available for special loads / applications. Contact LANCER for an Accessory Brochure that lists additional Rack Options.

Options

Depending on washer model, LANCER offers many options for your specific application including:

- Additional Dosing Pump(s)
- Drain Discharge Cool Down
- Steam Condenser
- Printers
- Steam Heating
- Seismic Restraints

Customization

Do you have a custom application? We offer customizable options including:

- Custom baskets and racks
- Custom chamber sizes (on certain models)
- Custom programming (pH Neutralization, etc)

Additional Services

LANCER offers additional services including:

- IQ/OQ Documents and Execution
- Preventative Maintenance Plans
- Extended Warranties
- User and Service Training

Cleaning Chemicals

Available exclusively from LANCER, our line of premixed detergents and acid neutralizers are specifically made for use in all high pressure washers. Direct dosing from bottle to washer via peristaltic pump ensures accurate cleaning.

LANCERCLEAN Detergent - 3 Formulas

- Potassium Hydroxide Based
- Sodium Hydroxide Based
- pH Neutral

LANCERACID Rinse - 3 Formulas

- Acetic Acid Based
- Citric Acid Based
- Phosphoric Acid Based



If you have any questions please contact LANCER.

Address:
3543 State Rd 419
Winter Springs, FL 32708
Website: www.lancer.com

Toll Free: (800)332-1855
Phone: (407)327-8488
Fax: (407)327-1229
Email: sales@lancer.com

JJD REV4 0408



Washer-Dryer Product Guide